sifin diagnostics gmbh | Berliner Allee 317-321 | 13088 Berlin | Germany | T: +49 30 700 144-0 | F: +49 30 700 144-30 | info@sifin.de | www.sifin.de

Anti-Salmonella O-Gruppen-Pools

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Anti-Salmonella O-Group Pools

Test Reagents for slide agglutination

INFORMATION FOR PROFESSIONAL USE

CE IVD

en

	O-Group-Pool	Vol.
Anti-Salmonella OMA	A, B, D, E, L	1 ml, 3 ml
Anti-Salmonella OMB	C, F, G, H	1 ml, 3 ml
Anti-Salmonella OMC	I, J, K, M, N, O, P	1 ml
Anti-Salmonella OMD	Q, R, S, T, U, V, W	1 ml

	O-Group-Pool	Vol.
Anti-Salmonella OME	X, Y, Z, 51, 52, 53	1 ml
Anti-Salmonella OMF	54-59	1 ml
Anti-Salmonella OMG	60-63, 65-67	1 ml

Intended use
The test reagents are intended for use in the serological identification of *Salmonella* strains isolated from test material of varied origin, using slide agglutination.
Additional serological differentiation must be performed using the group-specific test reagents and the Anti-Salmonella O, Vi test reagents in accordance with the White-Kauffmann-Le-Minor scheme (Kauffmann-White scheme).

Principle of the test
If the isolated strain possesses an antigen that is covered by the test reagent, this *Salmonella* antigen becomes bound when mixed with the specific antibody. The antigen-antibody reaction results in clearly visible agglutination of the strain.

Composition
The test reagents are mixtures of monoclonal antibodies. They are produced from cell culture supernatants of hybridoma cell lines which secrete antibodies against the respective group specific *Salmonella* O-antigens.
For the coverage of the O-Group K, the Anti-Salmonella OMC test reagent contains serum of immunized rabbits from which unspecific agglutinins have been removed by absorption.
Preservative: Sodium azide (NaN_3) 0.9 mg/ml)

Form in which product is supplied, shelf life and conditions for storage
The test reagents are available in liquid form ready to use.
If stored at 2...8 °C before and after opening, they may be used up to the date given on the label. After use, the bottle must be properly closed.
The test reagents may sometimes show turbidity which is not microbially caused. Such turbidity does not impair effectivity and the test reagents can be clarified by centrifugation or filtration.
The test reagents must have their temperatures adjusted to room temperature (18...26 °C) before use.

Warnings and precautions
The biotechnological manufacture of the test reagents means that the risk of contamination by infectious agents can be virtually excluded. Because they contain animal materials (fetal calf serum, stabilizer) and because of the addition of rabbit serum, they should be treated as potentially infectious and handled accordingly.
As these materials contain sodium azide, contact with the skin and mucous membranes must be avoided! In case of contact, rinse with plenty of water.
Since the performance of the slide agglutination test involves working with native pathogenic materials, all necessary work protection procedures must be adhered to (risk of infection)!

Materials and equipment not supplied
Glass slides, stirring rods, physiological saline, disposal containers for infectious material, White-Kauffmann-Le-Minor scheme.


Test material und methodology
The test material is streaked onto non-selective culture media (Nutrient Agar, Blood Agar) or Kligler Medium. Incubate for 16 – 20 hours at 35...37 °C.
Transfer a small amount of bacterial mass from a suspicious colony under test onto a slide and mix well with a drop of test reagent (ca. 25 µl), so that a homogenous, slightly milky suspension results. The slide should be placed on a dark surface.
The reaction is read with the naked eye by holding the slide in front of a light source against a black background and rocking it (tilting it back and forth).
To control for spontaneous agglutination, a negative control (NC) must be performed at the same time using physiological saline in place of the test reagent.

Test procedure
First carry out the test with the Anti-Salmonella OMA and Anti-Salmonella OMB test reagents. Approx. 98 % of *Salmonella* can be assigned to Anti-Salmonella OMA and Anti-Salmonella OMB.
If the strain does not agglutinate with these both test reagents, it is recommended to test this strain with Anti-Salmonella Vi, REF TR 1316. If this reaction is also negative the strain must then be agglutinated with the Anti-Salmonella OMC, Anti-Salmonella OMD, Anti-Salmonella OME, Anti-Salmonella OMF, and Anti-Salmonella OMG test reagents.

Evaluation
The test can only be evaluated if the negative control (NC) remains milky-opaque.
Positive: visible agglutination after the sample has been tilted back and forth less than 20 times. In a strongly positive reaction (PR), agglutination (coarsely or finely flocculent) appears as soon as the bacterial mass is mixed in. In a weakly positive result, agglutination only appears after the slide has been tilted back and forth 10-20 times.

Negative: if the suspension remains milky-opaque or the reaction begins to occur only after the slide has been tilted back and forth more than 20 times, the result is negative (NR).

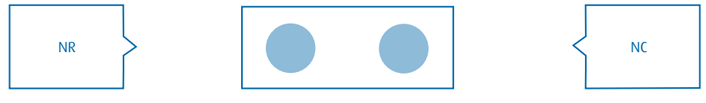Quality assurance during the testing procedure
For the quality control of serological determination of *Salmonella* type-O antigens by slide agglutination test, the good expression of the strain's cell surface antigen is important. The use of strains from interlaboratory tests, or field strains of defined origin that have been fully characterised by an external laboratory, or the use of suitable commercially available *Salmonella* test antigens is therefore recommended for quality control.

Limitations of the method
The test reagents react specifically with individual *Salmonella* strains, i. e. they only agglutinate *Salmonella* strains belonging to the O-Group Pools named in the declaration. In exceptional cases, there is the possibility of cross-reactions with other genera of *Enterobacteriaceae* (e. g. with strains of *Citrobacter* spp., *Hafnia alvei*, *Proteus* spp. or *E. coli*) due to antigen identities or related antigens. For this reason, biochemical examination is required to establish membership of the genus *Salmonella*.

Explanation of the symbols used

LOT	Batch code (Lot)
REF	Catalogue number
IVD	In Vitro Diagnostic Medical Device
TR	Test reagent
	Use by YYYY-MM (MM = end of month)
	Temperature limitation
	Consult instructions for use
	Slide agglutination

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