sifin diagnostics gmbh | Berliner Allee 317-321 | 13088 Berlin | Germany | T: +49 30 927 030-0 | F: +49 30 927 030-30 | info@sifin.de | www.sifin.de

Testseren Anti-Coli

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INFORMATION FOR PROFESSIONAL USE

en

Intended use
The test sera are used to test for the presence of typespecific antigens with *Escherichia* coli isolates from **piglets.** They are used as evidence of antigens to the cell surface of the *E. coli* isolate (O antigen, K antigen and F antigen) using slide agglutination (serotyping).

Principle of the test
If the *E. coli* strain possesses an antigen covered by the test serum, this antigen becomes bound when mixed with the specific antibody. The antigen-antibody reaction results in clearly visible agglutination of the strain.

Composition
The test sera are absorbed sera from immunised rabbits.
Preservative: Sodium azide (NaN_3) 0.9 mg/ml

Available products
Polyspecific: Anti-Coli P contains antibodies directed against all *E. coli* types listed as monospecific and the fimbrial antigen F4 (K 88).
Monospecific: Anti-Coli F 4 (K 88); Anti-Coli O 8:K 87; Anti-Coli O 138:K 81; Anti-Coli O 139:K 82; Anti-Coli O 141:K 85; Anti-Coli O 147:K 89; Anti-Coli O 149:K 91

Form in which product is supplied, shelf life and storage
The test reagents are provided in liquid form ready to use. If stored unopened at 2...8°C, the test sera may be used up to the date given on the label. After use, the bottle must be properly closed.
The test reagents may sometimes show turbidity not caused by microbes. Such turbidity does not impair effectivity and the sera can be clarified by centrifugation or filtration.
Prior to use, the test reagents must be warmed up to room tem-perature (18...26 °C).

Warnings and precautions
The test sera are biological materials (rabbit sera) which should be treated as potentially infectious and handled accordingly. As this product contains sodium azide, contact with the skin and mucous membranes must be avoided. In case of contact, rinse with plenty of water.
Since the performance of the slide agglutination test involves working with native pathogenic materials, all necessary work protection procedures must be adhered to (risk of infection)!

Materials and equipment not supplied
Distilled water, physiological saline (NaCl), glass slides, stirring rods, disposal containers for infectious material.


Test material and culture conditions
The test material is streaked onto non-selective culture media, such as Columbia Blood Agar, with the addition of sheep blood (50 ml /1 l) and is incubated for 16 - 20 hours at 35...37 °C. **Note: The F 4 (K 88) antigen is normally found in haemolytic types.**
The slide agglutination may only be performed with pure cultures (subcultures preferred).

Slide agglutination (procedure)
Transfer small amounts of bacterial mass from suspicious *E. coli* colonies onto a slide and mix each with one drop of the test serum (ca. 25 µl) so that a homogenous, slightly milky suspension results. Ensure that the slide is positioned on a dark surface.
The result is read with the naked eye by holding the slide in front of a light source against a black background and rocking it (tilting it back and forth).
To control for spontaneous agglutination, a negative control must be performed in parallel using 0.9 % NaCl solution in place of the test sera.

Interpretation of the results
The test can only be evaluated if the negative control (NC) remains milky-opaque.
Positive: visible agglutination after the sample has been tilted back and forth less than 20 times. In a strongly positive reaction (PR), agglutination (coarsely or finely flocculent) appears as soon as the bacterial mass is mixed in. In a weakly positive result, agglutination only appears after the slide has been tilted back and forth 10-20 times.

Negative: if the suspension remains milky-opaque or the reaction begins to occur only after the slide has been tilted back and forth more than 20 times, the result is negative (NR).

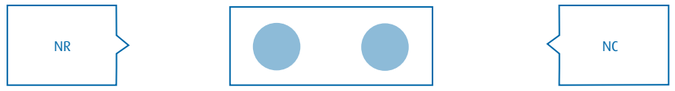Test procedure
First carry out a screening agglutination with the polyspecific test serum Anti-Coli P. If the reaction with the polyspecific test serum is positive, the isolate must then be typed with the monospecific test sera.
Evidence of fimbrial antigen: Anti-Coli F 4 (K 88);
Evidence of possible O-K types: Anti-Coli O 8:K 87; Anti-Coli O 138:K 81; Anti-Coli O 139:K 82; Anti-Coli O 141:K 85, Anti-Coli O 147:K 89, Anti-Coli O 149:K 91
Isolates that react with the polyspecific antiserum or with the monospecific antiserum Anti-Coli F 4 (K 88) but not with any of the monospecific O:K antisera can in principle by typed but display an antigen profile that deviates from the standard and which cannot therefore be verified.

Quality assurance during the testing procedure
For the quality control of serological determination by slide agglutination test, the good expression of the strain's cellsurface antigen is important. The use of strains from interlaboratory tests, or field strains of defined origin that have been characterised by an external laboratory, is therefore recommended for quality control.

Limitations of the method
The test reagents react with *E. coli* strains, which contain antigens of the specificity named in the declaration.
In exceptional cases, there is the possibility of crossreactions with other *E. coli* antigens or other *Enterobacteriaceae* due to related antigens.

Explanation of the symbols used

LOT Batch code (Lot)
REF Catalogue number
Use by YYYY-MM (MM = end of month)
Temperature limitation
Consult instructions for use

Date of revision: 01/11/2017

sifin diagnostics gmbh | Berliner Allee 317-321 | 13088 Berlin | Germany | www.sifin.de